UPS Stool Collection Instructions

Preparation


1. Before collecting your sample, put on the gloves included in the kit. Write your name and date of birth on both tube labels and the biohazard specimen bag. Write your date of birth below your name on the tube with the liquid.

2. Open the paper collection bowl and follow the instructions to install it on your toilet. Have a bowel movement into the bowl.

Sample Collection

3. Grab the tube **without** liquid and use the scoop attached to the lid to put stool into the tube, filling it between halfway and the fill line. Next, remove the tube **with** liquid from the small plastic bag, discard the bag, and use the scoop attached to this lid to add stool into this tube until the stool or the liquid reaches the fill line. **Do not overfill the tubes. Screw the caps on both tubes tightly.**

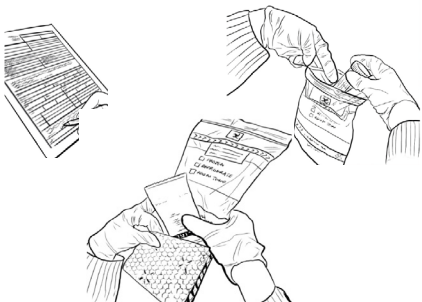4. Write the date and time of collection on **BOTH** tubes, the biohazard bag, and the Diarrhea Pathogen Requisition or test order/paperwork received from your provider. Place **BOTH** tubes (not the paperwork) in the biohazard bag and seal it. If the ice pack is not ready, refrigerate the samples. **DO NOT place samples in the freezer.**

5. Place **BOTH** the ice pack and the biohazard bag in the silver pack. Seal the silver pack.

Send Sample via UPS

6. Place the sealed silver pack and the paperwork in the pre-paid return shipping pack. Seal the shipping pack. Document the tracking number located on the label of the pre-paid shipping pack by either taking a picture of the label or writing the tracking number on the UPS Patient Quick Reference Guide.

7. Drop off the shipping pack at a UPS store or drop box, or schedule a pickup via UPS. You can visit www.ups.com or call 1-800-742-5877 to find a store or drop box location or schedule a pickup. If you schedule a pickup, please reference the following information:

Account Number: 7A280X Account Number Zip Code: 07755

Other

Adhesive-remover pad

8. If needed, use the adhesive remover pad to remove adhesive from the toilet bowl.

Please do not send anything other than what has been provided to you in this kit.

If you have any questions, please feel free to contact Genesis Laboratory via phone at 732-389-1530, option 0, or email at genesisresults@genesislaboratory.com

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